

# inset® 30 Quick Reference Guide

inset® 30




inset® 30 combines insertion device and infusion set into a single sterile, portable unit. Simply follow the steps below.




Check the tamper-proof seal and sterile paper—do not use device if the seal or paper has been broken. Pull the red tape to remove the seal.



Pull off the sterile paper. Do not touch the spring.

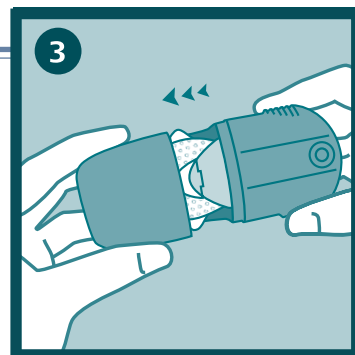


Hold on the lined indentations with one hand and pull the lid off with the other.\*  
**CAUTION:** Be careful not to bend or touch the inset® 30 introducer needle.

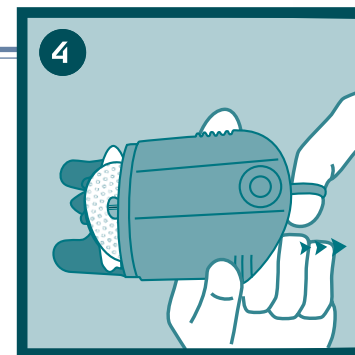

Place your fingers on the lined indentations on each side. Pull back the spring softly until you hear a “click”. The needle protector will lift up or peel off. Do not touch or bend the needle.




Immediately before insertion, carefully bend back the needle protector and remove it. This will remove part of the adhesive backing paper. Be sure the tape is not stuck on the introducer needle.

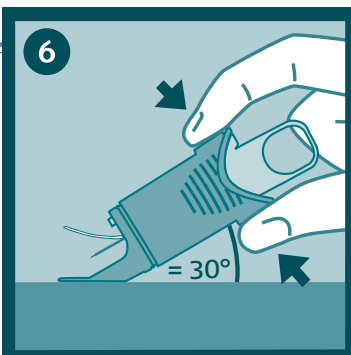


**INSERTION:** Place your index finger on the top release button and your thumb on the bottom release button. Position inset® 30 with its legs level (flat) on the skin to ensure an insertion angle of 30°. Make sure to maintain a 30° angle while inserting and gently squeeze once to insert inset® 30.



Make sure that the cannula stays in place: Put your finger gently on the clear window while removing the insertion device and introducer needle. Pull it gently straight back.

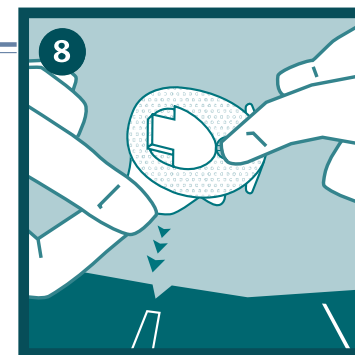


Secure the cannula with one finger on the clear window. Carefully remove the backing paper under the cannula housing as shown.

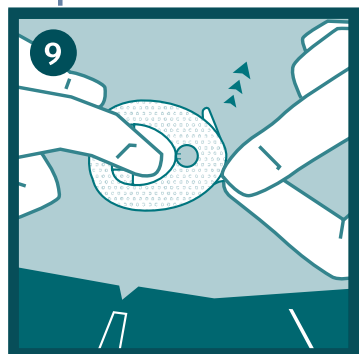

Remove the backing paper under the front end. Massage the adhesive thoroughly on to the skin.  
**IMPORTANT:** Do not attempt to straighten the adhesive if curled.



Put the lid back on for disposal.

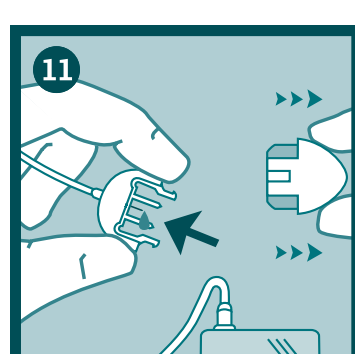

Remove the circular protective cap. Prime infusion set until medication dispenses from the tip. Ensure there are no air bubbles in the tubing or inset® 30. Consult your pump manufacturer’s instructions for use for proper priming procedures. Place a finger in front of the cannula housing and insert the connector needle. Make sure the connector locks with a “click.” Fill the empty space in the cannula with the following prime:

U100 insulin:  
0.7 units (0.007mL)

**\*IMPORTANT NOTE!** If the infusion set is not securely placed in the insertion device with the needle pointing straight ahead, pain or minor injury may occur during insertion.

Please note that this guide only illustrates some of the steps to follow when inserting inset® 30 and it is not a substitute for reading the instructions for use which comes with the inset® 30 product.

When initially trying a new infusion set, we recommend this to be done under the supervision of a healthcare provider.

inset® is a registered trademark of Unomedical a/s.

Manufactured by:

